

# The HARVEY

## Stethoscope

**Two Kinds of Physicians  
Will Want the Harvey  
Elite Cardiology  
Stethoscope. . .**

**Those Who Need the Best.  
Those Who Want the Best.**

### The Harvey tradition

The Harvey stethoscope was invented by Dr. W. Proctor Harvey, one of the world's leading cardiologists and experts on auscultation. The Harvey name is synonymous with quality and acoustical effectiveness. Now, Welch Allyn/Tycos offers a new Harvey cardiology stethoscope that is so comfortable and durable, it will quickly become the most essential tool in your busy practice.




### Exceptional acoustic effectiveness

Most significant heart sounds occur in the 200-500 Hz frequency range, while human auditory sensitivity is limited below 1000 Hz. The new Harvey Elite Stethoscope provides superior acoustical performance in this most important range. Its quality stainless steel chest-piece offers a bell for transmission of low-frequency sounds (20- 100 Hz) and a diaphragm for detection of higher-frequency heart and lung sounds (100-1000 Hz).

In addition, the Elite's dual-bore, one-piece PVC tubing provides separate sound pathways to each ear for unrivaled accuracy in sound transmission, without loose tubes that might create ambient noise. Whatever your specialty, you'll be able to diagnose difficult heart and lung sounds more accurately and effectively than ever before.

### Unexcelled user comfort and portability

The Harvey Elite Stethoscope is so lightweight and comfortable to use and carry, you'll find it an indispensable diagnostic companion, even over extended periods. Its quality stainless steel binaurals, hand-crafted to meet Welch Allyn/Tycos' exacting specifications, are fully rotatable to fit the individual user's ear canal. Each Elite comes with two sets of soft and hard ear tips, threaded for secure connection and easy removal, to ensure user comfort and eliminate disturbing environmental noise. The Elite also features a non-chill rim on both the bell and diaphragm to allow comfortable application to your patient's skin.

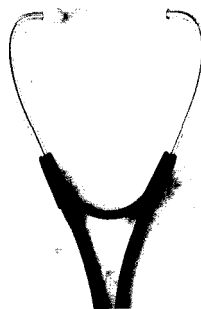

### Durable construction throughout

The Harvey Elite's flexible PVC tubing lasts longer than latex and folds to fit easily in the practitioner's pocket. All metal parts are precision-crafted of stainless steel for years of trouble-free service. This fine instrument also features a long-lasting, triple-leaf binaural spring that can be adjusted for proper tension and comfortable fit through years of use.

### Pediatric Capability

The Harvey Elite Stethoscope offers you the flexibility of pediatric capability in the same instrument you use on older patients. Standard accessories include both pediatric bell and diaphragm chestpieces for effective auscultation of neonates and infants - at no additional cost.

### Warranty

The Harvey Elite is backed by Welch Allyn/Tycos' unparalleled customer service. Our 10 year warranty on parts and construction augments the respected Harvey tradition in auscultation with Welch Allyn/Tycos' industrywide reputation for quality.

### Ordering Information

All Welch Allyn/Tycos stethoscopes and accessory parts can be purchased from an authorized distributor in your area. Call Welch Allyn directly for the name of the dealer nearest you.

### Complete Elite Stethoscopes

5079-122 22" Black 5079-270 25" Burg.  
5079-125 25" Black 5079-271 25" Navy

### Elite Accessories/Replacement Items

#### Eartip Pairs

Soft	Hard
5079-110 Large	5079-113 X-Large
5079-111 Medium	5079-66 Large
5079-233 Comfort Sealing	

#### Tubing

5079-119 22" Black 5079-272 25" Burg.  
5079-120 25" Black 5079-273 25" Navy

#### Binaural/Spring

5079-121 Black 5079-274 Burgundy  
5079-275 Navy

#### Diaphragm Assembly

Pediatric	Adult
Black 5079-116	Black 5079-278
Gray 5079-279	Gray

#### Bell Non-chill Ring

Pediatric	Adult
5079-117 Black 5079-118	Black 5079-276
Gray 5079-277	Gray

#### Diaphragm Non-chill Rim

Pediatric	Adult
5079-127 Black 5079-126	Black 5079-185
Gray 5079-184	Gray

#### Diaphragm Disc

5079-124 Pediatric 5079-123 Adult

#### Harvey Elite

5079-232 Accessory Kit (Black)  
5079-280 Accessory Kit (Gray)

Black Burgundy Navy

Welch Allyn/Tycos Inc.  
95 Old Shoals Road  
Arden, NC 28704  
Tel: (704) 684-4895  
Fax: (704) 687-1002

**Welch Allyn**