


SALTER LABS®



*Three-Channel
Oxygen Supply Tubing*
*Superior safety due to crush-resistant internal
wall design...a Salter innovation*

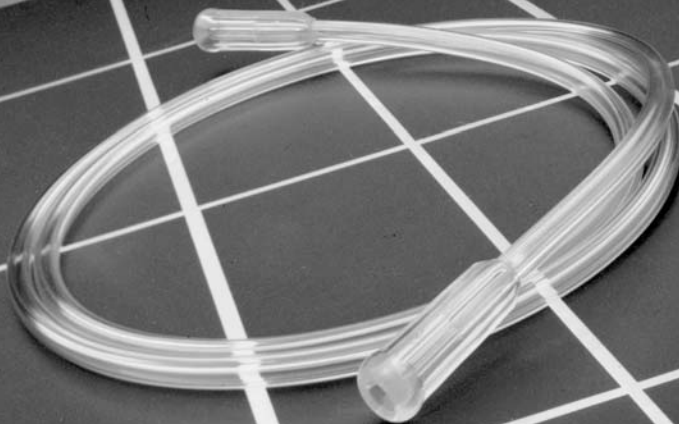

Features

- Three-channel safety tubing
- When uncoiled, tubing remains straight
- No “memory”, clear tubing
- Innovative connector design.
“Ribbed” end piece with rounded
and tapered edges.

Benefits

- Oxygen supply cannot be occluded or interrupted
- Virtually eliminates “torque” or pull on cannula
- Will not coil-up on floor. More aesthetically pleasing
- Easy attachment to all fittings. Helps to eliminate
entanglement with furniture or door jambs. No square edges



2002 Three-Channel Oxygen Supply Tubing

Special design has three inside channels that help prevent the tubing from becoming occluded. Ideal for home patients where long lengths can cause a potential kinking problem. Available in lengths from 4' to 100'. *Special or custom lengths available... a Salter service oriented benefit.*

Note: Oxygen supply tubing also available in smooth bore configuration.

Innovative Tubing "End Fitting"*

- Vinyl material with sufficient clarity to view patency
- Designed with "Ribbed Body" to permit easy gripping when making connections
- Rounded and tapered edges to prevent snagging or catching on various surfaces

NEW, redesigned end fittings "ribbed" for sure grip




rounded front with tapered rear edges

three channel safety tubing



End Fitting

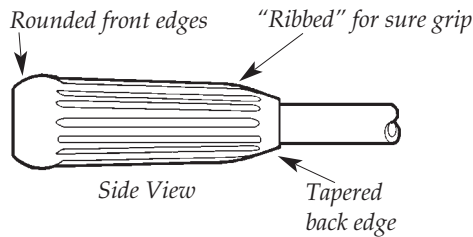

Ordering Information

Oxygen Supply Tubing - 3/16" I.D.

	Units (per case)	Safety Channel Clear	Smooth Bore Clear	Safety Channel Green
4' oxygen supply tube	50	2004	2504	
7' oxygen supply tube	50	2002	2503	2002G
9' oxygen supply tube	50	2009	2509	2009G
10' oxygen supply tube	50	2010	2510	
14' oxygen supply tube	50	2014	2514	2014G
15' oxygen supply tube	50	2015	2515	2015G
20' oxygen supply tube	25	2020	2520	2020G
21' oxygen supply tube	25	2021	2521	2021G
25' oxygen supply tube	25	2025	2525	2025G
30' oxygen supply tube	20	2030	2530	
35' oxygen supply tube	20	2035	2535	2035G
40' oxygen supply tube	20	2040	2540	2040G
50' oxygen supply tube	20	2050	2550	2050G
75' oxygen supply tube	15	2075	2575	
100' oxygen supply tube	10	2100	2500	

Special lengths available on request. Minimum order is one (1) case.

Safety Oxygen Supply Tubing with Thread Grips

7' safety oxygen supply tube with one molded connector end, one female thread grip end	50	2099	
7' safety oxygen supply tube with female thread grip on both ends	50	9997	
7' safety oxygen supply tube with male thread grip on both ends	50	9998	
7' safety oxygen supply tube with one male thread grip end, one female thread grip end	50	9999	
50' safety oxygen supply tube with one molded connector end, one female thread grip end	20	2099-50	

Special lengths available on request. Minimum order is one (1) case.

Special Tubing

1/8" I.D. clear tubing (200' coil)	1	1018
3/16" I.D. white tubing (100' coil)	1	1019

SALTER LABS®

100 W. Sycamore Road - Arvin, CA 93203
 Phone: 661-854-3166
 1-800-421-0024 (CA only)
 1-800-235-4203 (U.S. & Canada)
 Fax: 661-854-3850
 Toll Free Fax: 1-800-628-4690 (U.S. & Canada)
 International Fax: 661-854-3850
 E-mail: salterlabs@us.salterlabs.com
 www.salterlabs.com

CE

0482

ISO 9001/EN 46001 & ISO 13485